The Benchmark for Reliability

## **Compass** Ambidextrous Treatment Centres



### **Compass** delivery options :

Unrivalled console mobility and Compass's determined versatility bring all the benefits of convenience, flexibility, efficiency and Patient comfort.

Belmont's base-mounted delivery system design is an understated advancement in dentistry and Patient care.

Ambidextrous capability across the range doubles the qualification of Compass in any configuration to perform excellent service for any Operator and procedure.

COMPASS 'Surgery System' (Ambidextrous) Air or Electric version

COMPASS 'Flexible Dr System' (Ambidextrous) chair and Operator's console




COMPASS 'Flexible System' (Ambidextrous) chair and cuspidor with Assistant's console



### for better daily dentistry :

From an elegant, slim base-mounted arm comes a responsiveness which ensures Operators achieve comfortable working positions.

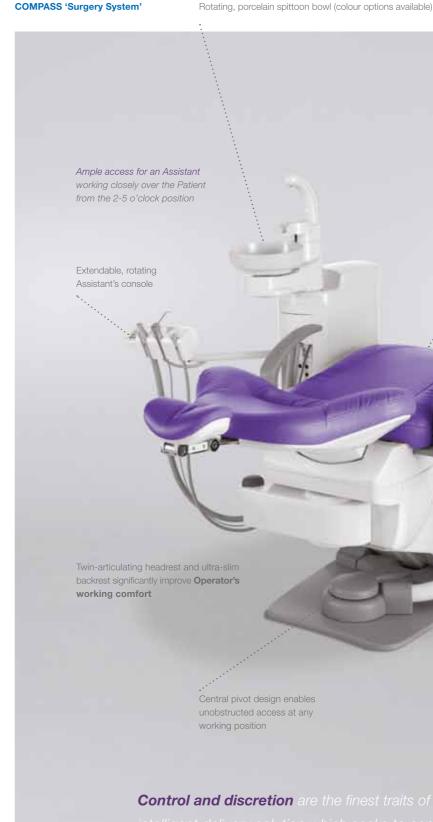



Luxury upholstery (Special order item)

Operator console key features ~ 4 programmable Operator settings LED digital display First-priority handpiece selection Digital timer 'E' version additional features ~ Preset micro-motor Torque control with Auto-Reverse 2 fibre-optic Turbine hoses

> Operator's console rotates smoothly on a base-mounted arm designed for precise and balanced positioning

Control and discretion are the finest traits of the Compass delivery system. An

### Increased surgery efficiency & Patient access :



Base-mounted swing-arm delivery system provides the full range of Operator access

To improve efficiency and increase workflow are the key aims of the Compass treatment centre. Sheer mobility enables simultaneous tasks to be carried out, thus saving the Patient and surgery time.



## comfort for all Dentists & Patients :

Correct ergonomics and smooth movements ensure that Compass is a pleasure to work with at any position.

COMPASS 'Surgery System'

Extensive access for Assistants working at 2-5 o'clock

# Simulataneous preparation and consultation saves valuable working time between every Patient. The uniquely perfect clean and preparation arrangement to the rear of the chair also improves lines of communication between Operator and Assistant. Clean & Prep zone Consultation zon or detach spittoon bowl

60 second changeover from left- to right-side delivery, or vice versa requires no tools or expertise. The simple procedure results in identical conditions for left- and right-handed Operators.






Right-handed set up. Raise chair height and remove the rear cover.

Incline chair. 1-2. Swing table round to left side.



Compass

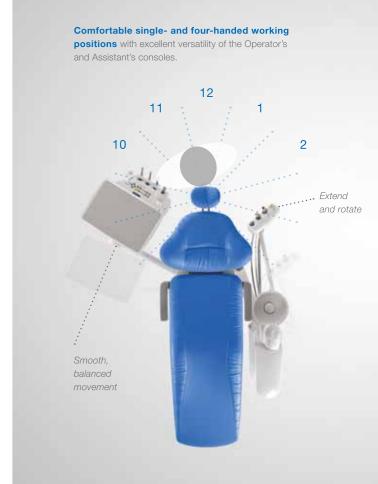



3-4. Swing cuspidor round to right side.



Rear cover is easily removed/replaced without tools or expertise



5. Secure cuspidor in postion. Secure cuspidor locking lever. 2. Rotate console. Left-handed set up.

## Compass 'Flexible' treatment centres :

*Flexible configurations to suit any Practice and procedure. The flexibility of Electric (E) and Air (A) versions to fulfil individual preferences.* 

### for your dentistry :

Tailor the treatment centre for your exact preferences.



Operator's control panel includes instrument, chair, spittoon and light controls



'E' version controls additionally include preset micro-motor to run torque control with Auto-Reverse, and 2 fibre-optic Turbine hoses

Assistant's console with chair, spittoon and light controls





Ambidextrous cuspidor with rotating and detachable porcelain spittoon bowl

<sup>1</sup> Light controls on 'E' version only



ON/OFF Touchless sensors on all lights








700 SERIES LIGHT

Chair- or ceiling-mounted<sup>1.</sup> **4,200 Kelvin, Halogen light** Three preset intensities: 18,000 Lux, 28,000 Lux and Composite Curing Safe Mode

Soft-push button with intensity indicator features on 700 Series and 820S models

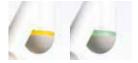

820S LIGHT Chair-mounted 5,000 Kelvin, Halogen light Three preset intensities: 18,000 Lux, 28,000 Lux and Composite Curing Safe Mode

Fluid, triple-axis head movement features on 820S and 900 Series models

### 900 SERIES LIGHT

Chair- or ceiling-mounted **5,000 Kelvin, ten LEDs** Stepless intensity control: 4,000 Lux to 32,000 Lux Composite Curing Safe Mode (5,500 Lux)



LED bulbs consume less power and maintain significantly lower heat during operation

Spittoon bowl colour options: blue/green/pink glass

Control all chair movements by hand or foot: by hand from both the Operator's and Assistant's console or by using the foot-control. Spittoon functions and ON/OFF light control <sup>1</sup> are also activated from both the Operator's and Assistant's console.





Foot control for instruments and chair

# Belmont

### takara.co.uk/dental

### Takara Belmont (UK) Ltd

London head office & showroom Belmont House, One St. Andrews Way, Bow, London E3 3PA UNITED KINGDOM Telephone +44 (0) 20 7515 0333 Facsimile +44 (0) 20 7987 3596 Email dental@takara.co.uk

Manchester showroom Telephone +44 (0) 161 743 9992

Showroom visits by appointment







The manufacturer reserves the right to alter specifications without prior notice.